



FreeBrowse: User Guide v1

User Interface

Landing Page: Click on “Get started” to create a new project

Welcome to
Freebrowse

Get started

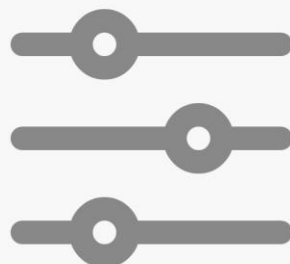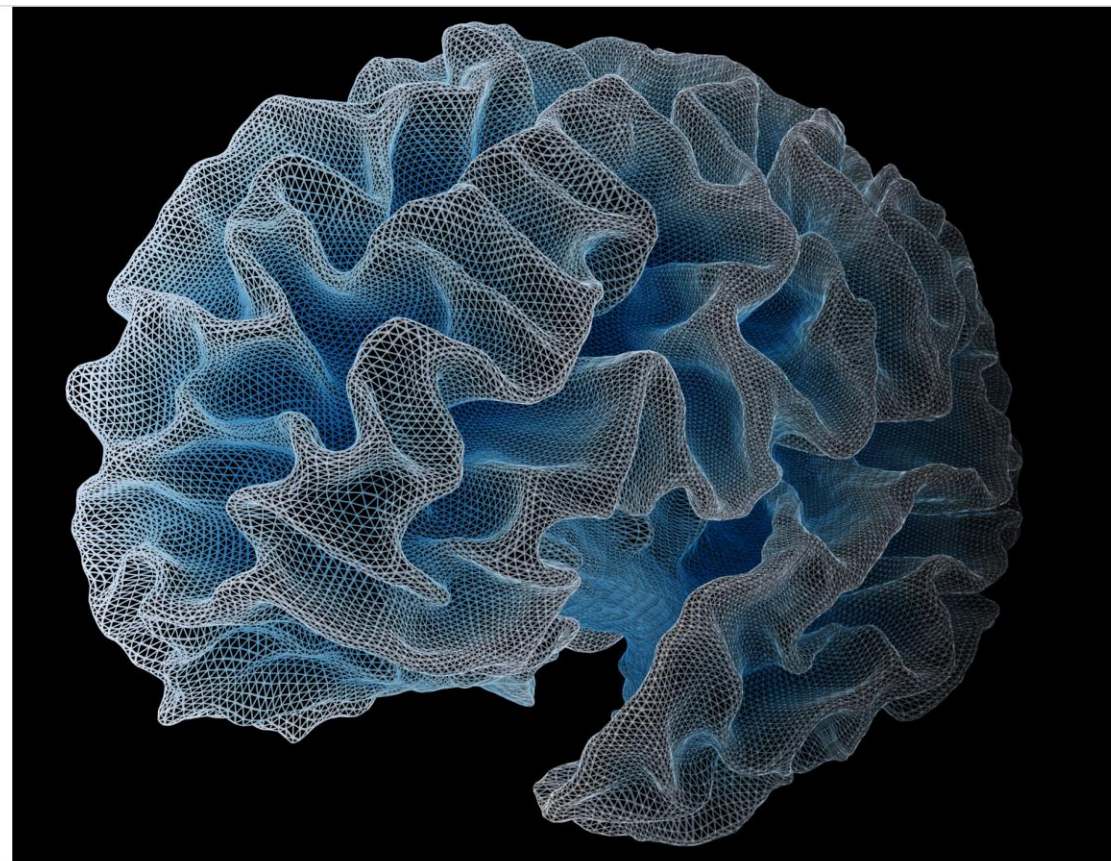


FreeBrowse is a browser based version of FreeSurfer's Freeview


This is an early proof-of-concept prototype. It is under development and will change in the upcoming weeks and months. Please send feedback (bugs, feature requests and anything else) to pwighton@mgh.harvard.edu

Creating a new project:

Choose Project Name and select volumes and surfaces for loading



Load volumes & surfaces




My computer

Cloud

Project Name:

Subject_1



Drop files to attach, or [browse](#)

Cancel

Open

FreeBrowse Interface – Areas Explained


Toolbar - Viewports - Left Column - Right Column

Toolbar: Main interactions in FreeBrowse
(change mode, change layout, share /
download / save project, access user profile, etc.)

Left Column:
Details about the
loaded data
(Loaded files,
file settings &
voxel info)


Viewports:
Interactive viewer
(sagittal, coronal,
axial, 3D Brain)

Right Column:
Manipulating the
loaded data
(Brush settings,
Commenting)

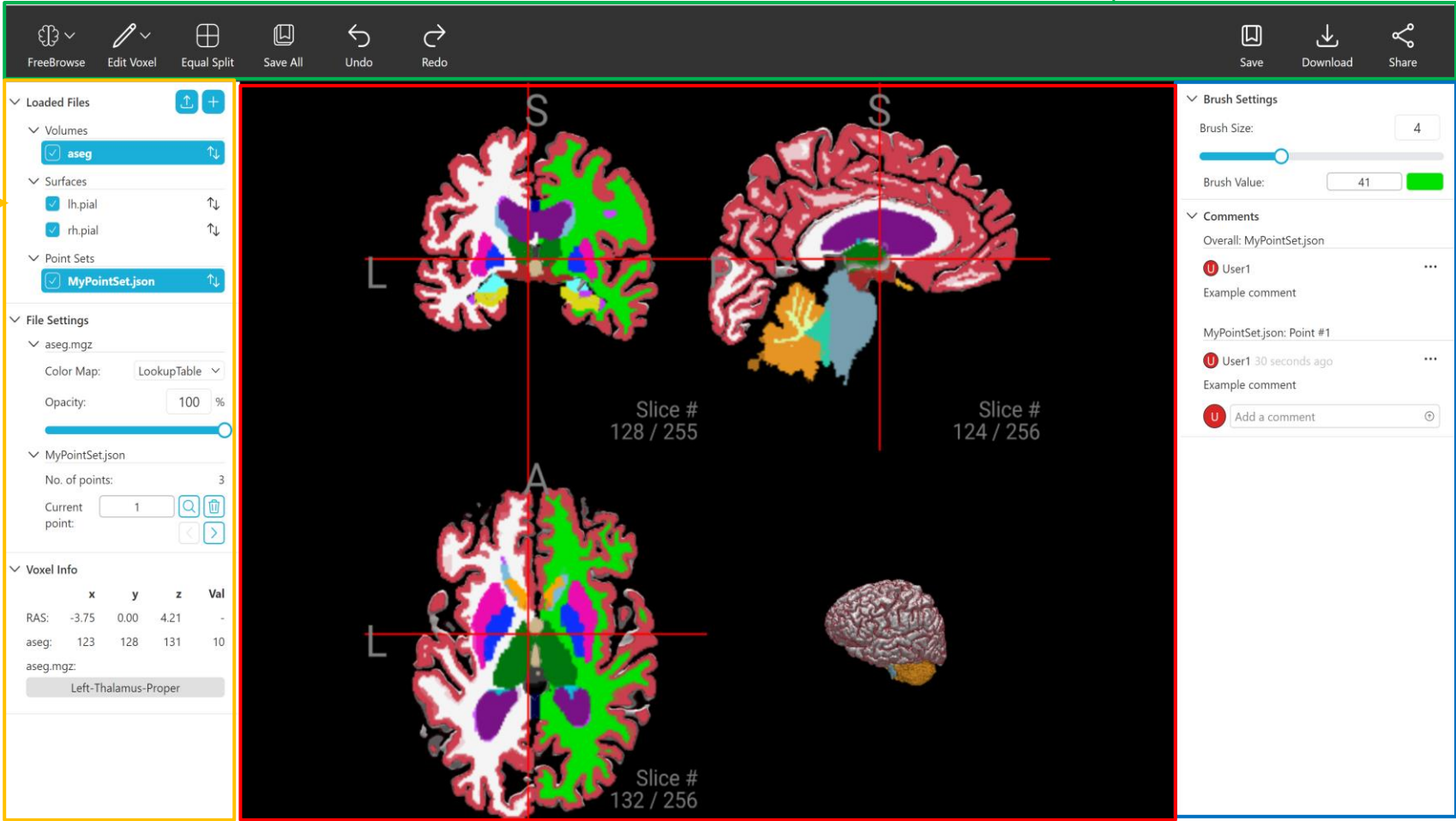

Navigate viewports

- Hover over the view you want to select
- Scroll through the slices – Ctrl + scroll or Up and Down arrow keys
- Zoom in & out – Scroll

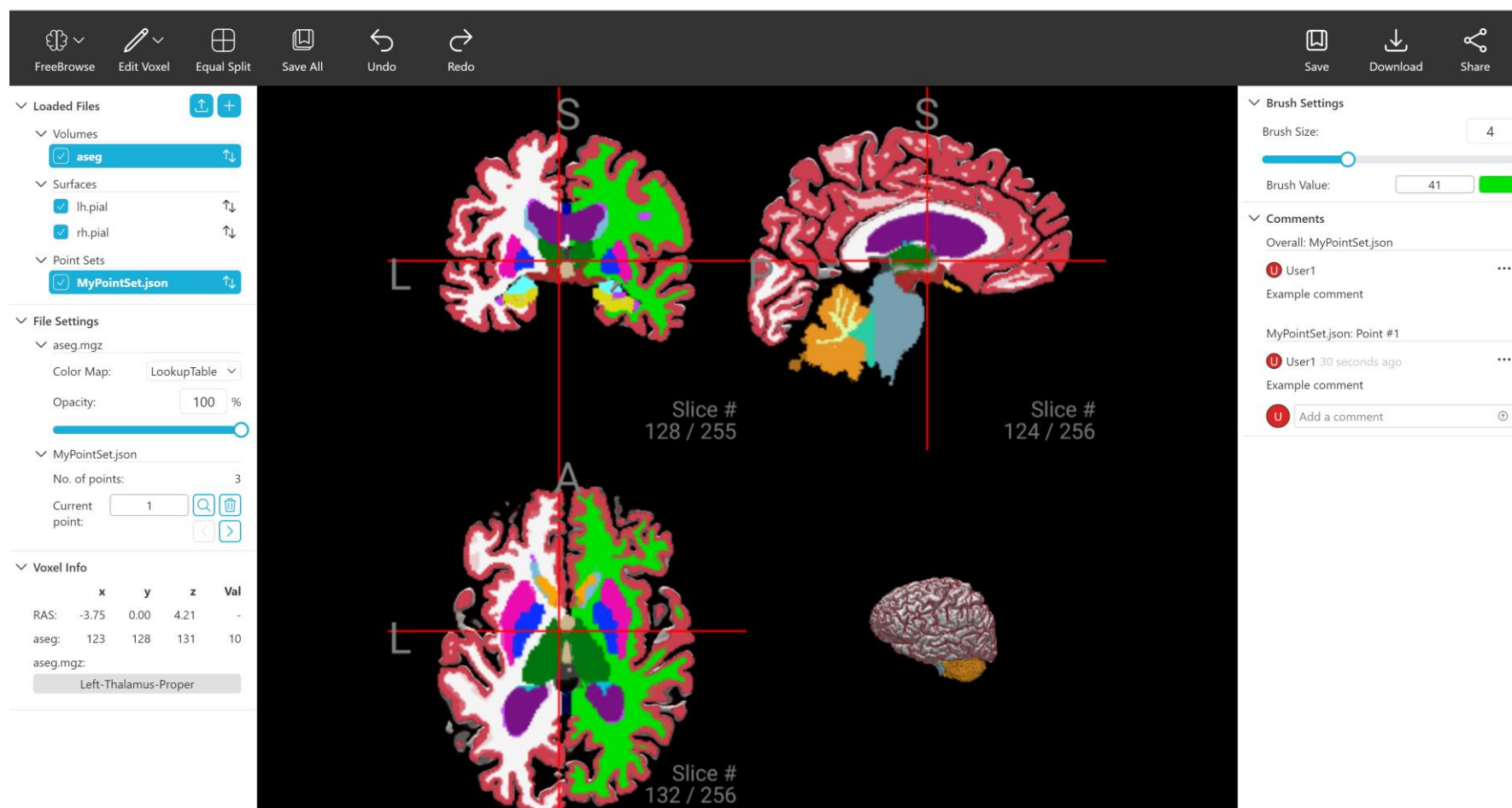

Adjust file settings

- Select the file that you want to interact with and use the controls in the File Settings area

FreeBrowse

Edit Voxel

Equal Split

Save All

Undo

Redo

Save

Download

Share

Loaded Files

Volumes

aseg

Surfaces

lh.pial

rh.pial

Point Sets

MyPointSet.json

File Settings

aseg.mgz

Color Map: LookupTable

Opacity: 100 %

MyPointSet.json

No. of points: 3

Current point: 1

Voxel Info

	x	y	z	Val
RAS:	-3.75	0.00	4.21	-
aseg:	123	128	131	10
aseg.mgz:	Left-Thalamus-Proper			

Brush Settings

Brush Size: 4

Brush Value: 41

Comments

Overall: MyPointSet.json

User1

Example comment

MyPointSet.json: Point #1

User1 30 seconds ago

Example comment

Add a comment

View voxel information

- Click on the desired voxel in the viewport
- Find information like the coordinates, the assigned label or the voxel value in the left column

FreeBrowse

Edit Voxel

Equal Split

Save All

Undo

Redo

Save

Download

Share

Loaded Files

Volumes

aseg

Surfaces

lh.pial

rh.pial

Point Sets

MyPointSet.json

File Settings

aseg.mgz

Color Map: LookupTable

Opacity: 100 %

MyPointSet.json

No. of points: 3

Current point: 1

Voxel Info

	x	y	z	Val
RAS:	-3.75	0.00	4.21	-
aseg:	123	128	131	10
aseg.mgz:	Left-Thalamus-Proper			

S

L

A

Slice # 128 / 255

Slice # 124 / 256

Slice # 132 / 256

Brush Settings

Brush Size: 4

Brush Value: 41

Comments

Overall: MyPointSet.json

User1

Example comment

MyPointSet.json: Point #1


User1 30 seconds ago

Example comment

Add a comment


Create a new point set

- Click on the “+” button next to Loaded files
- Enter a name for your point set & choose a color

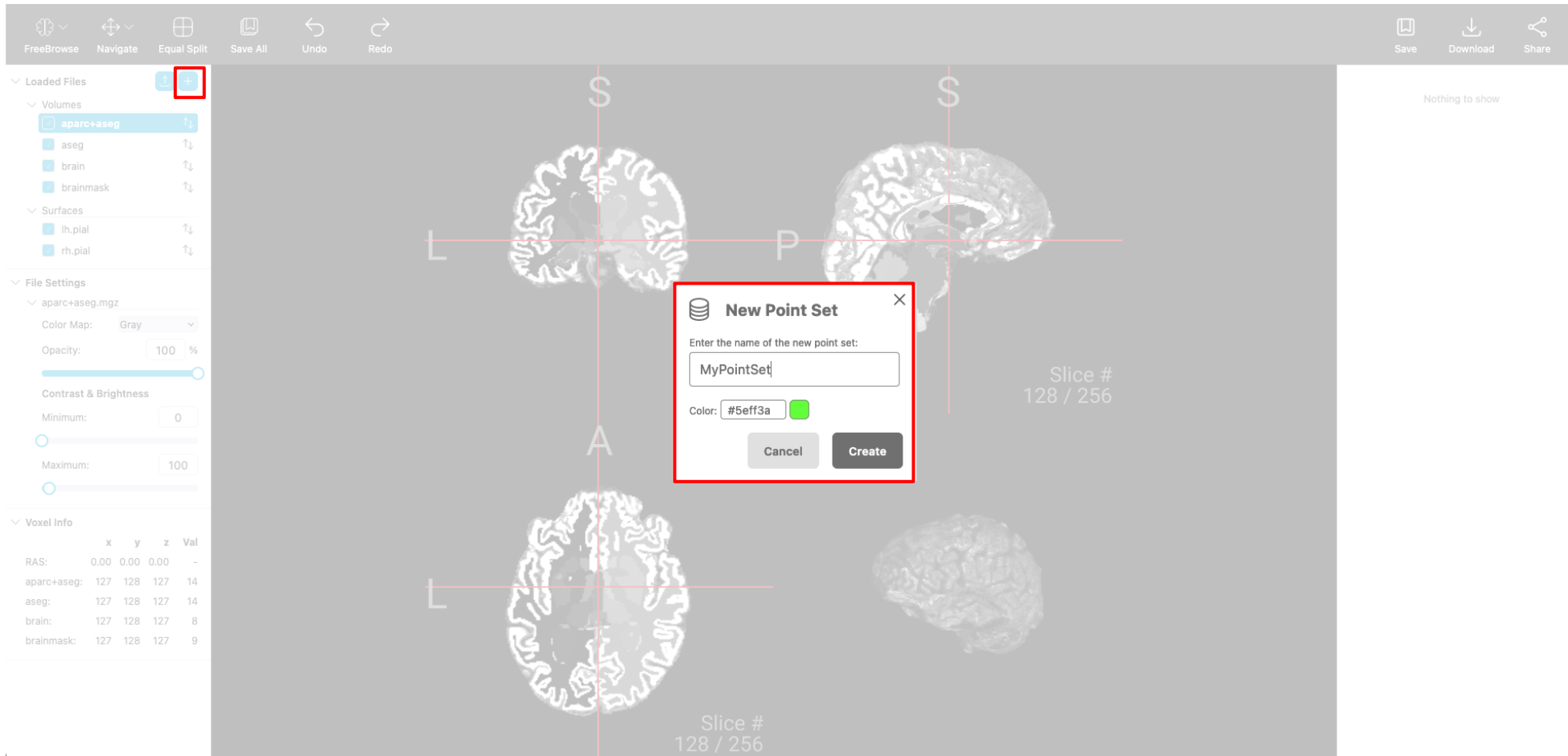


Create & navigate between waypoints

- Select the “Edit Points” mode
- Select the point set that you want to interact with
- Create new points with left click
- Navigate between points using the controls in the “File Settings” area

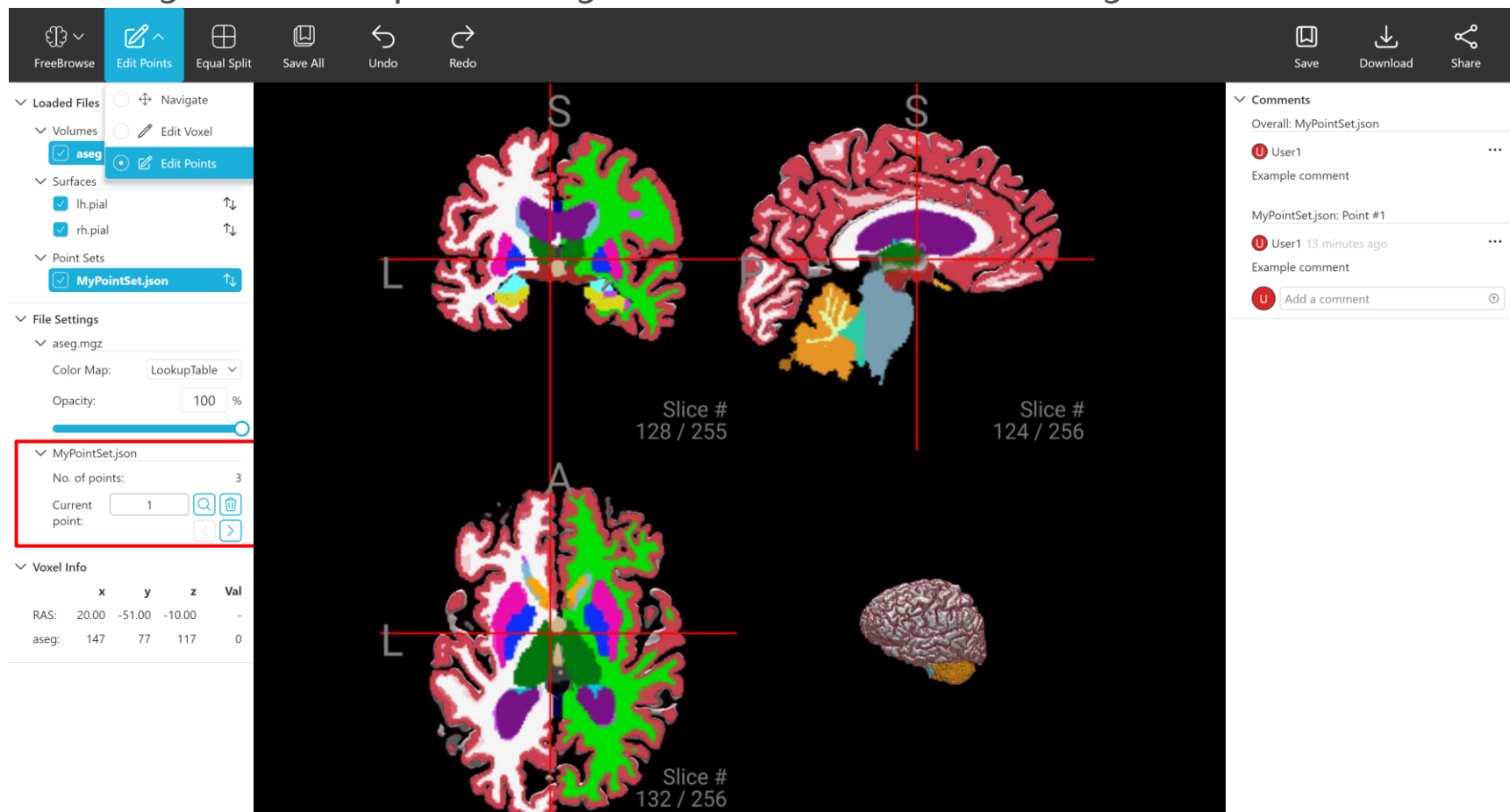

Edit Voxels: Change the brush settings, draw over the slices to edit the voxels

- Select the “Edit Voxel” mode
- Set the preferred brush size and value from the “Brush Settings” area
- Press and hold the left mouse button to edit voxels
- When you are done press the “Save” button

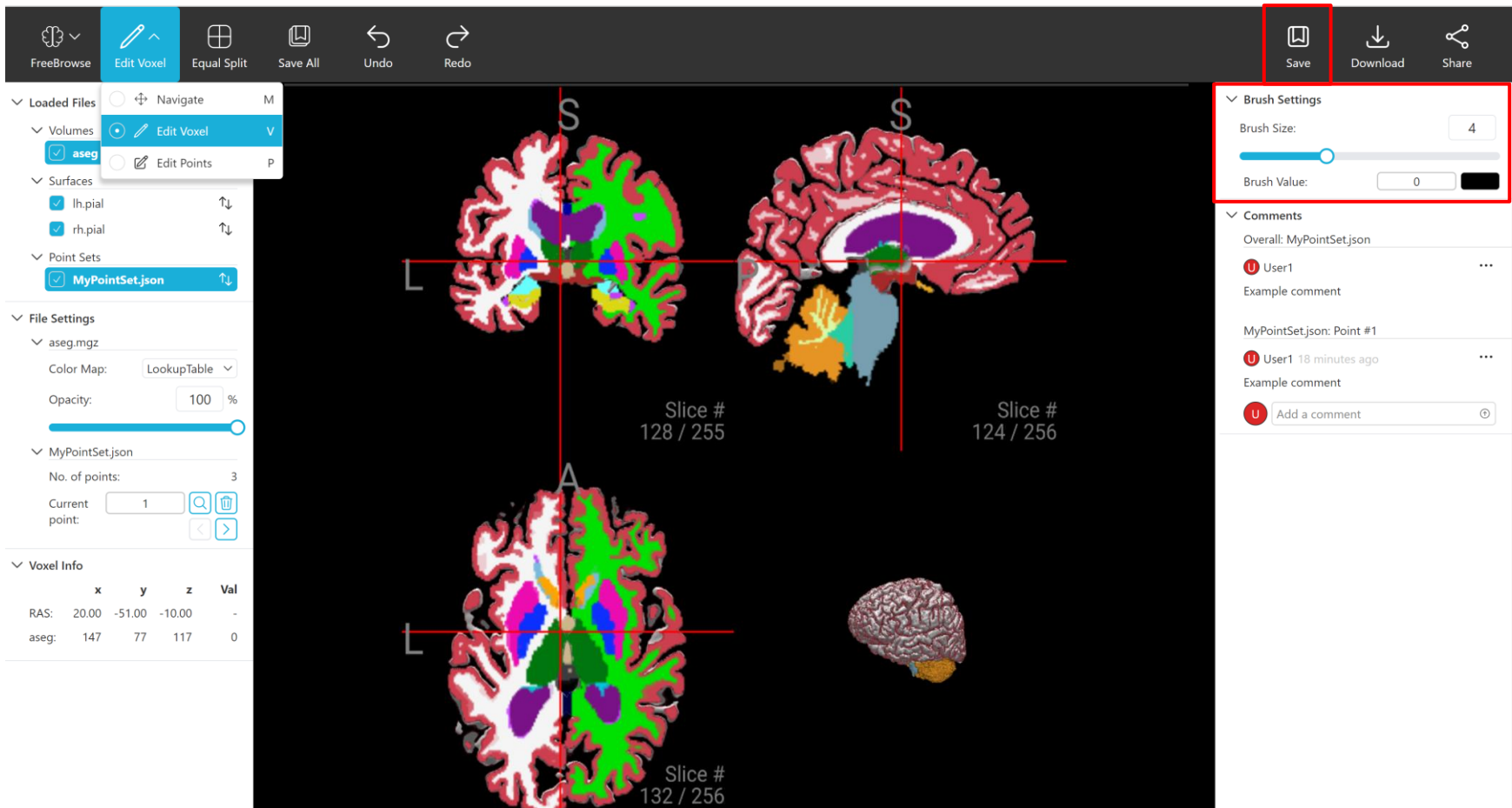

Add comments to waypoints

- Navigate to the waypoint that you want to comment
- From the “Comments” area you can comment on
 - The entire pointset
 - A single waypoint

FreeBrowse

Edit Voxel

Equal Split

Save All

Undo

Redo

Save

Download

Share

Loaded Files

Volumes

aseg

Surfaces

lh.pial

rh.pial

Point Sets

MyPointSet.json

File Settings

aseg.mgz

Color Map: LookupTable

Opacity: 100 %

MyPointSet.json

No. of points: 3

Current point: 1

Voxel Info

	x	y	z	Val
RAS:	-3.75	0.00	4.21	-
aseg:	123	128	131	10
aseg.mgz:	Left-Thalamus-Proper			

S

L

A

Slice # 128 / 255

Slice # 124 / 256

Slice # 132 / 256

Brush Settings

Brush Size: 4

Brush Value: 41

Comments

Overall: MyPointSet.json

User1

Example comment

MyPointSet.json: Point #1


User1 30 seconds ago

Example comment

Add a comment

Deeplink: Share the exact same view with another person for collaboration

- Navigate to the exact view that you want to share
- Press the “Share” button
- This will copy a link to your clipboard that you can share

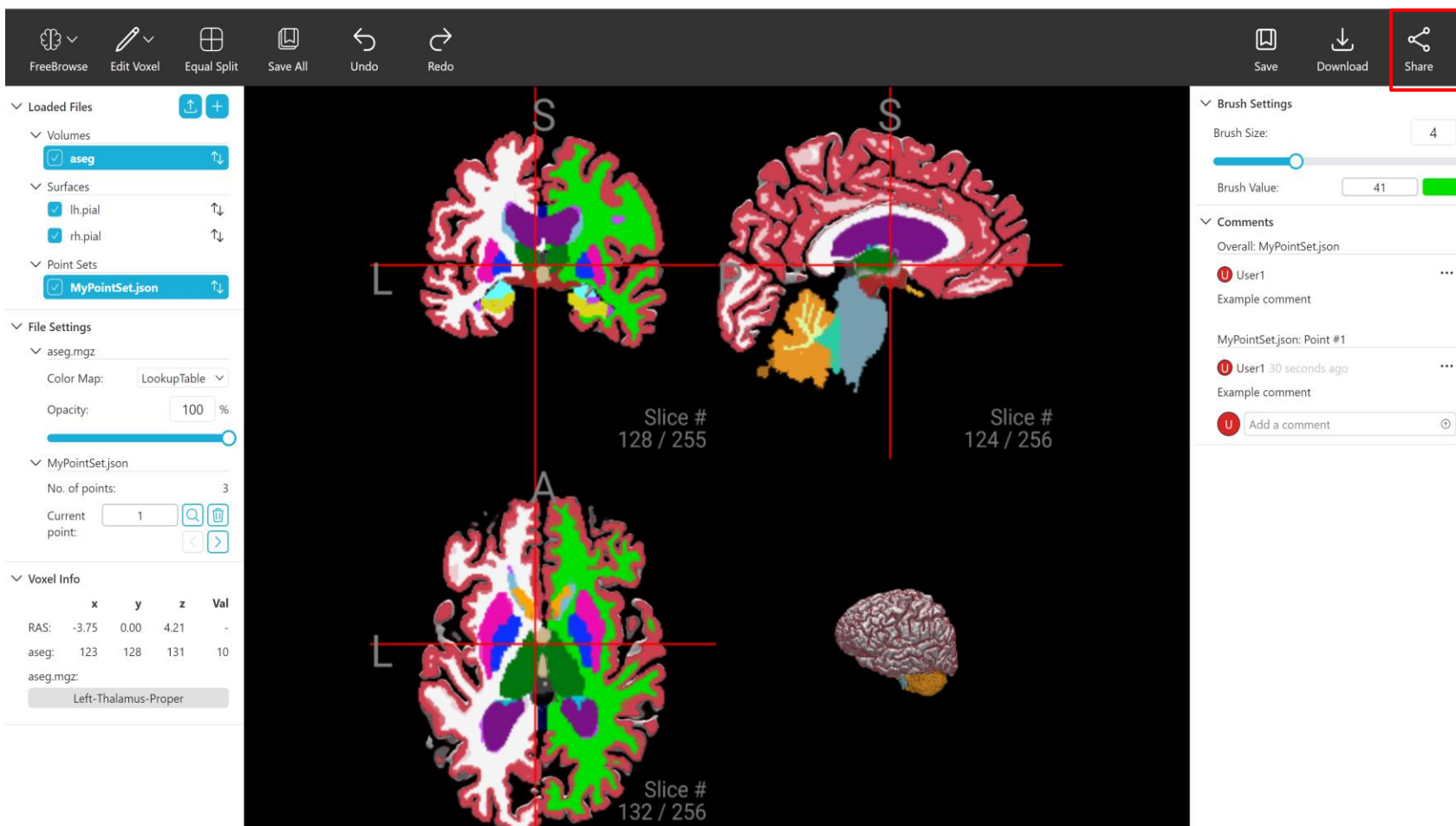

Summary of Interactions

Interactions & Shortcuts

FreeBrowse is a **single page application** that can be controlled with a **mouse** or **trackpad** and a **keyboard**.

PoC Interactions:

Desired action	FreeBrowse interaction	Shortcut (Windows)	Shortcut (MAC)
Select	Left-click with mouse	-	-
Zoom	Scroll with mouse	Scroll	Scroll
Rotate	Left-click hold & drag mouse across 3D viewport	-	-
Panning	Right-click & drag mouse across 2D viewport	-	-
Undo	Left-click “Undo” button	CTRL + Z	CMD + Z
Redo	Left-click “Redo” button	CTRL + Y	CMD + Y
Navigate through slices	Arrow Up / Down	CTRL+ scroll with mouse/ up and down arrow keys	CMD + scroll with mouse/ up and down arrow keys
Save	Left-click “Save” button	CTRL + S	CMD + S
Draw	Left-click hold & drag mouse across 2D viewport (“Voxel Edit” mode only)	-	-
Create waypoint	Left-click with mouse (“Point Edit” mode only)	-	-